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Venus Pro II Hair Removal Device

[Important Taboos]

- Pregnant women, individuals with skin diseases, those sensitive to sunlight, and those wearing medical electronic devices should avoid use
- Not suitable for use on damaged, swollen, or tattooed areas
- Do not look directly at the light source; avoid shining it into your eyes
- Not waterproof. Keep away from water sources and humid environments

[Hair removal Management Plan]

- Weeks 1-4: Use 3-4 times per week
- Weeks 5-8: Use 2-3 times per week
- After Week 8: Use once every half month or once a month

* Please use it appropriately based on the growth of hair. For heavy hair, you can increase the gear and frequency accordingly. The specific usage varies from person to person.

[Gear and mode]

- Multiple gear settings available, starting from the lowest gear when first using
- Refer to the table below to determine which light intensity is appropriate for your skin and hair color
- Supports multi-area adaptation for face, armpit, arms, legs, bikini, etc.
- Equipped with functions for ice application, skin rejuvenation, and hair removal, allowing for switching between modes

Gear usage reference table (skin color card)

Skin Color						
Hair Color						
Dark Golden / Light Brown	1-6	1-6	1-5	1-4	×	×
Brown	1-6	1-6	1-5	1-5	×	×
Dark Brown	1-6	1-6	1-5	1-5	×	×
Black	1-6	1-6	1-6	1-5	×	×

(Note: The skin color marked with the symbol "X" in the above table is not suitable for using this product.)


[Quick steps]

Use a razor to remove the hair from the surface of the skin and clean it. (It is recommended to use lubricating foam in conjunction.)

Connect the power supply and turn on the device. Rotate the dial from 1 to 6 to start the operation. The ice point mode is automatically enabled.

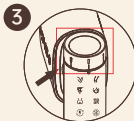Select the corresponding skin care mode and energy level.

Before using the hair removal device, please wear protective goggles to protect your eyes from harm.

Align the light outlet with the skin and press the light button to start hair removal. (You can also enable the continuous flashing mode for sliding hair removal.)

After the hair removal is completed, you can choose the SR Photon Skin Rejuvenation mode to apply another treatment, which will promote the skin's smoothness and texture, completing the skin care process.